


RANGE OF ACCESSORIES



SPECIFICATIONS				
Modelo	Ref. 010098	Ref. 010099	Ref. 010103	Ref. 010107
Height (mm)	150	150	170	160
Width (mm)	180	170	170	148
Depth (mm)	180	180	180	180
Weight	0,540 kg	0,689 kg	0,735 kg	0,697 kg

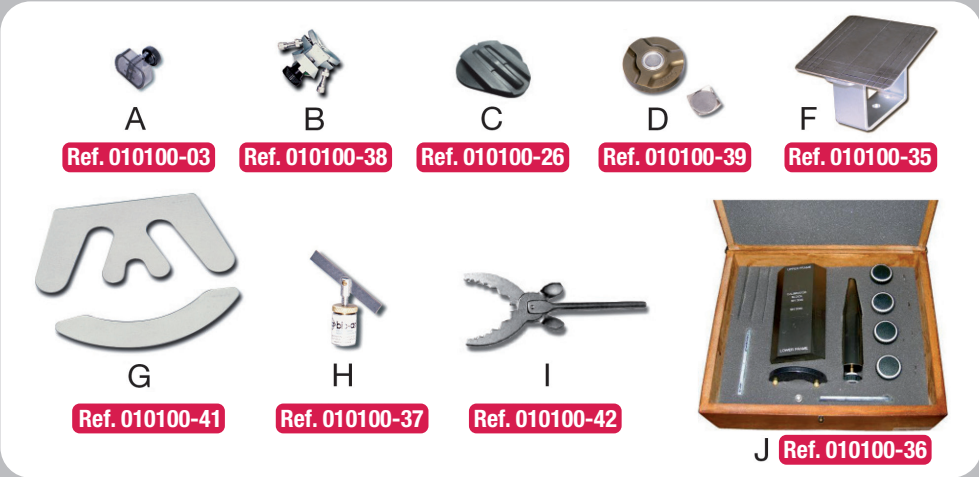

Ref. 010099

Articulator + Case



Ref. 010103

Articulator + Case



Ref. 010107

Articulator + Case



Ref. 010098

Articulator + Carboard box

The 4 models are perfectly compatible with the Standard and Élite type facebows.

ORIGINAL ARTICULATORS

Ref. 010098 - 010099 - 010100-55 - 010102-50 - 010103 - 010107





## Arcon type articulator “ARQUÍMEDES PROFESIONAL” (semiadjustable) A7 PLUS



Robust structure made of aluminum, ensuring greater stability. Design that provides greater visibility.

- ✓ Semiadjustable / Arcon.
- ✓ Fixed intercondylar distance (110 mm).
- ✓ Curve condylar guide.
- ✓ Regulable condylar guide.
- ✓ Regulable Bennet angle.
- ✓ Central lock.
- ✓ Stabilization system of the condilar guide movement with silicone connection.
- ✓ Support for upper branch in open position.

Ref. 010099



bio-art  
intelligent solutions

## “ARQUÍMEDES JUNIOR” (semiadjustable) A7 FIX

- ✓ Semi-adjustable / Chest.
- ✓ Fixed intercondylar distance (110 mm).
- ✓ Condylar guide fixed at mid 30°
- ✓ Bennet angle fixed at mid 15°
- ✓ Curved condylar guide.
- ✓ Magnetic stabilization system (embedded) for movement of the condylar guide.
- ✓ Pin for upper ramus support in open position.

Ref. 010098

bio-art  
intelligent solutions



## Arcon type articulator “ARQUÍMEDES PROTRUSAL” (semiadjustable) A7 Plus-E



Similar to the Archimedes Pro model, it has exclusive aluminum condylar guides with millimeter adjustment.

- ✓ Semiadjustable / Arcon.
- ✓ Fixed intercondylar distance (110 mm).
- ✓ Regulable condylar guide.
- ✓ Regulable Bennet angle.
- ✓ Flat condylar guide.
- ✓ Bilateral locking system.
- ✓ Support for upper branch in open position.
- ✓ Stabilization system of the condilar guide movement with spring.

Ref. 010103




bio-art  
intelligent solutions

## Standard face-bow

The Standard face-bow is an apparatus used to record the spatial position of the patient's upper dental arch in relation to the skull and transfer the record to the articulator.

Transfer of the patient record is performed by setting the complete face-bow in the upper frame of the articulator.

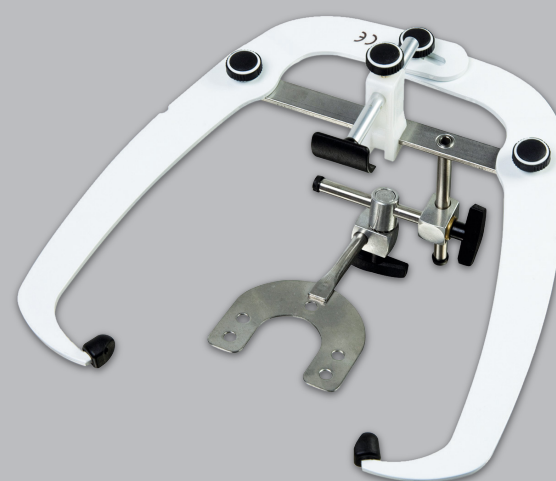

Ref. 010102-50

bio-art  
intelligent solutions

CE Medical Device

## Arcon type articulator “4.000 S” (semiadjustable) 4.000 S



- ✓ Semiadjustable / Arcon.
- ✓ Adjustable intercondylar distance (expanding spindle).
- ✓ Curve condylar guide.
- ✓ Regulable condylar guide.
- ✓ Regulable Bennet angle.
- ✓ Stabilization system of the condilar guide movement with O-ring.

Ref. 010107



bio-art  
intelligent solutions

## Élite face-bow

The new Élite Face-Bow brings an articulated and fixation system for the fork that allows a faster registration process for the professional and comfortable for the patient. It has a jig transfer assembly that eliminates the needs for face-bow in the assembling process of the casts in the articulator, turning the process faster, more practical and stable.

Its Naziun Relator with adjustable height, including an infra-orbital indicator, turns the registration process more precise. The face-bow is equipped with a sliding system which allows smooth lateral movement (opening/closing) and facilitates the overall use. The ear pieces are anatomical, removable and chemically sterilized.



Ref. 010100-55

bio-art  
intelligent solutions

CE Medical Device