



# PONI CHAIRS

Refs. 100600 - 100602



# PONI CHAIR without backrest

Ref. 100600

The chair incorporates a safety braking system. When no one is sitting, the wheels lock to prevent the chair from moving accidentally. When someone sits down, they unlock and rotate freely.

## Specifications:


Maximum seat height: 730 mm

Minimum seat height: 530 mm

Depth: 320 mm • Width: 360 mm

Weight: 7 kg • Maximum load: 150 kg

## Wheels with safety brakes (EN 1335)

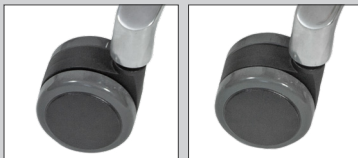

# PONI CHAIR with backrest

Ref. 100602


## MULTIPLE ADJUSTABLE ROTATION

The backrest can be swayed slightly independently.

The whole backrest can be raised and lowered and free.


The air rod can be raised and lowered and rotated.

The pulley can rotate 360°. Soft tone, no damage to the floor.

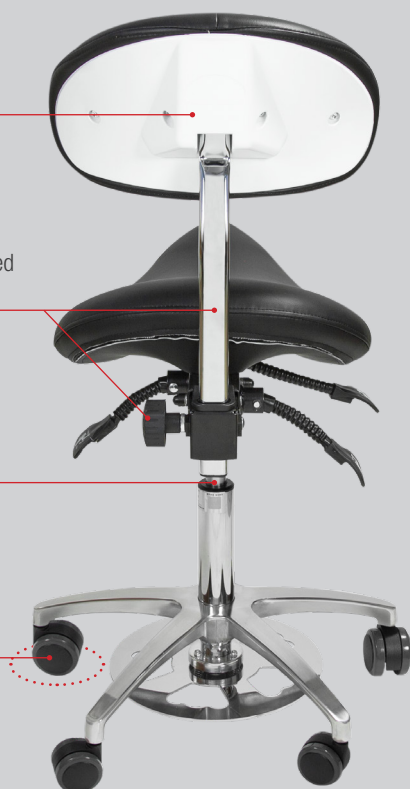

## ROCKING FUNCTION

## Wheels with safety brakes (EN 1335)

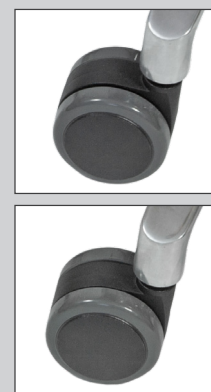

You can control the lift with your foot.



- ✓ **Height Adjustable.**
- ✓ **High quality PU leather.**
- ✓ **Effective against disinfection of treatment rooms with everyday cleaning agents.** Not easy to fade and age.
- ✓ **High resilience foam.** Rigorous control of raw materials. High than the national standard of high resilience, to bring a better sitting feeling.

- ✓ **Four-stage inlet pneumatic bar.** Safe + Smooth.TUV and SGS Authentication.
- ✓ **Hydraulic oil seal.** Multiple protection, no oil leakage, no air leakage! Long service life.
- ✓ **Oil filter cleanliness.** Built-in oil scraping device. The protruding part is oil-free and clean, without fear of false contact.

## Specifications:

Maximum seat height: 640 mm

Minimum seat height: 200 mm

Depth: 320 mm • Width: 370 mm

Weight: 11 kg • Maximum load: 150 kg



**MESTRA**

**Talleres Mestraitua SL**  
(+34) 944 530 388  
info@mestra.es - www.mestra.es  
48150 SONDIKA (Spain)



Management System  
ISO 9001:2015

www.tuv.com  
ID 1100018003