



MESTRA

CAPSULE MIXER

Ref. 100625

In modern dental practice, the durability and success of a restoration critically depend on the homogeneity of the initial mix. Manual dosing or unstable vibration seriously compromise the mechanical properties of the material.

The **capsule mixer** has been designed under strict standards to provide a high-speed oscillatory motion -controlled by a microprocessor- which guarantees a porosity-free mix with the perfect gloss and exact setting time. In short: a material with professional quality, optimal for its use.

The most significant features of this device are:



HIGH SPEED

Adjustable vibration frequencies of up to 5000 rpm, ensuring a homogeneous mix and perfect consistency.



UNIVERSAL COMPATIBILITY

Designed for capsules of different materials. A single piece of equipment for multiple restoration solutions.



PRECISION DIGITAL CONTROL

Equipped with an LCD screen and tactile membrane buttons for quick timer and speed configuration. Furthermore, cleaning is extremely simple.



ERGONOMIC AND SAFE DESIGN

Its advanced internal spring damping system drastically reduces vibrations transferred to the workbench while maintaining a low noise level. This improves both patient and professional comfort.




STABILITY AND QUIET OPERATION

Its advanced internal spring damping system drastically reduces the vibrations transferred to the table of work while maintaining a low noise level. This improves the comfort of the patient and the professional.



DURABILITY

Manufactured with a robust housing and heavy-duty components for a long service life.

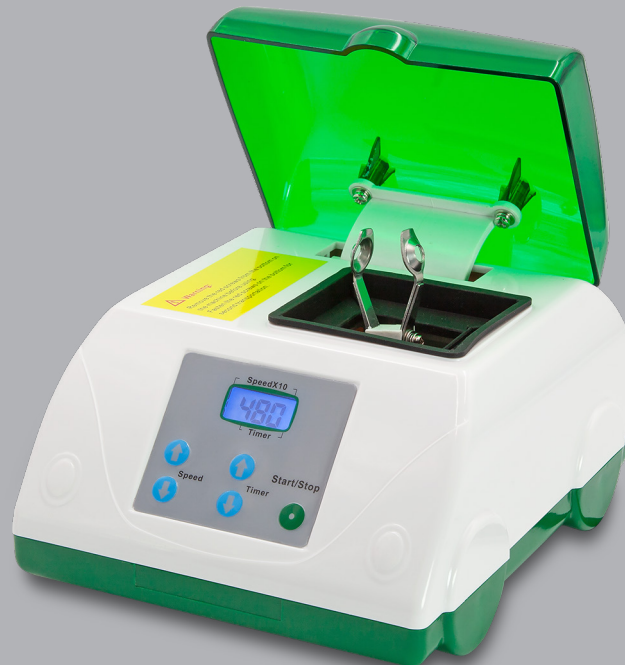

Characteristics:

Dimensions: 210 x 175 x 150 mm	Timer: 3 - 180 "
Packaging dimensions: 255 x 215 x 185 mm	Fuse: 2 A
Net weight: 1,87 kg	Power: 40 W
Gross weight: 2,1 kg	Tension: 220 V, 50 Hz
Speed: 2800 - 5000 rpm	Noise: ≤ 65 dB (A)

PACKAGE CONTENT:

- ✓ 1 mixer.
- ✓ 2 replacement fuses.
- ✓ 2 replacement springs.



MESTRA

Talleres Mestraitua SL
(+34) 944 530 388
info@mestra.es - www.mestra.es
48150 SONDIKA (Spain)



Management
System
ISO 9001:2015

www.tuv.com
ID 1100018003